



## THE ROHO<sup>®</sup> AirLITE<sup>®</sup> CUSHION

### COMFORT

### DURABILITY

### POSITIONING

### LIGHTWEIGHT

### SIMPLICITY

The ROHO AirLITE is unlike any cushion you have ever seen. Designed for comfort and peace of mind, the AirLITE combines the benefits of a sealed air support pad and a durable contoured foam to provide postural support and optimal durability. The non-adjustable pre-set air support pad, uses patented ROHO Air Flotation technology, while a contoured high density foam base positions the body in a comfortable position, creating one of the most effective and easy-to-use cushions on the market today. And at only 2 pounds, it continues the tradition of lightweight products from The ROHO Group.

AirLITE...comfort, durability and positioning in one simple package...the solution that others can't match.



# AirLITE® CUSHION

## COMFORT

By combining a contoured, anatomical-shaped foam with a pre-set, sealed Air Flotation component, the AirLITE enables clients to sit comfortably for an extended time. AirLITE's technology and design provides the comfort clients want and the support they need.

## DURABILITY


Foam wears out. Most foam cushions will compress over time. The AirLITE cushion's pre-set, sealed Air Flotation component provides comfort and peace of mind you can count on long after other "just foam" cushions have lost their effectiveness.

## POSITIONING

The contoured shape of the AirLITE cushion provides clients with support for basic positioning needs. Leg (femur) troughs and an ischial well provide support for lower extremity control, stability and comfort.

## LIGHTWEIGHT

The AirLITE cushion weighs only 2.17 pounds (1kg). Compare this to alternative technologies at two to four times the weight of the AirLITE. A lightweight cushion makes a mobility system easier to propel and transport.

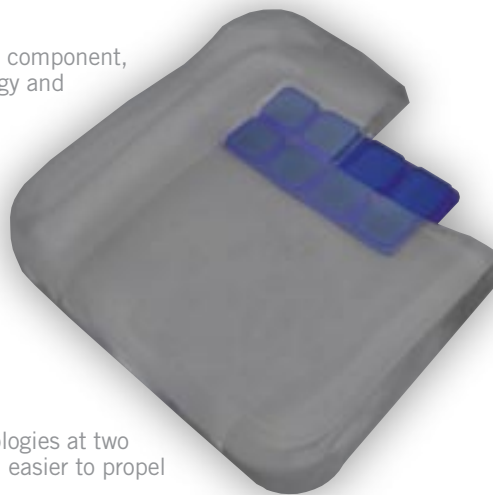

## APPLICATIONS

### THE AirLITE CUSHION IS APPROPRIATE FOR CLIENTS WHO:

- have partial to full sensation
- are able to perform independent weight shifts
- require basic posture support and comfort
- require basic positioning needs to assist in the prevention of skin breakdown
- need a lightweight, non-adjustable seating system

## SPECIFICATIONS

### INCLUDES:

AirLITE Cushion, Instruction Booklet, cover

### CONSTRUCTION:

Contoured Foam base – Flame Resistant Polyurethane Foam

Pre-Set, Sealed Air Cell: Polyurethane

Cover: Fluid resistant polyurethane film laminated to 2 way stretch Polyester knit, and non skid bottom.

### SIZES:

Available to fit chairs: 15" x 15" / 38cm x 38cm, 16" x 16" / 41cm x 41cm, 16" x 18" / 41cm x 46cm, 17" x 17" / 43cm x 43cm, 18" x 16" / 46cm x 41cm, 18" x 18" / 46cm x 46cm, 20" x 18" / 51cm x 46cm

### WEIGHT:

Varies by size, approximately 2.17 lbs. / 1 kg (18" x 16")

### WEIGHT LIMIT:

250 lbs. / 114 kg

### WARRANTY:

24-month limited

### MEDICARE PART B CODE:

E2605, Positioning Wheelchair Seat Cushion

## FEATURES & BENEFITS

### SIMPLICITY:

It is easy. With the pre-set Air Flotation component, just put the AirLITE cushion in the wheelchair and comfort, enhanced alignment and durability are achieved. No adjustment necessary.

### FLUID RESISTANT:

The breathable, fluid resistant cover protects the cushion from odor and is machine washable.

### SECURITY:

The non-skid bottom increases security of the cushion in the chair.

The AirLITE is designed for clients with tactile sensation and who are low risk for pressure ulcers. For clients who are at high risk or have existing pressure ulcers, we recommend the use of ROHO DRY FLOATION® Cushions.

## MORE PRODUCTS FROM ROHO



CONTOUR SELECT™ CUSHION




PRODIGY MATTRESS OVERLAY®



MOSAIC® CUSHION

### FOR ADDITIONAL INFORMATION, PLEASE CONTACT:

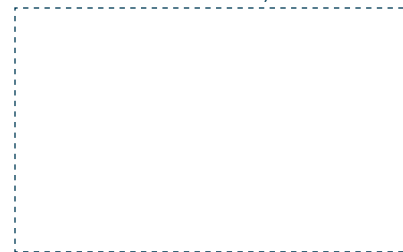

100 North Florida Avenue • Belleville, IL 62221-5429 • USA • [www.therohogroup.com](http://www.therohogroup.com)  
U.S. & Canada: 1 800 851 3449 • 1 618 277 9173 • fax 1 888 551 3449  
Outside the U.S. & Canada: 1 618 277 9150 • fax 1 618 277 6518

©2007 ROHO, Inc. The following are registered trademarks of ROHO, Inc.: ROHO®, shape fitting technology®, DRY FLOATION®, AirLITE®, CONTOUR SELECT®, PRODIGY Mattress Overlay®, and MOSAIC®.

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend specifications presented in this sell sheet.

LITG02006  
SCHW 8M  
3/08